Gasoutlet Panel. Product Description.



Product Description	The Gasoutlet Panels are suitable in preparation rooms, emergency rooms and wherever the space is limi- ted on the walls. To save space the terminal units are mounted diagonaly abowe and below each other. The panels are available for ressesed and exposed installation. The terminal units are self closing ie. closing automaticaly when connected equipment is disconnected. The terminal units are according to hospital standard, this means that equipment can not be connected to wrong terminal unit. The gas outlets are furnished with a service valve which means that the gas flow will be closed when the outlets are removed for service. The terminal units in the panel are delivered with gaslabels according to standard and protection plug. The Gasoutlet Panel is equiped with built in evacuation ejektor for anesthesia gases, the capacity is adjustable up to 25 lpm.		
Recycling	The product shall be recycled according to local regulations.		
Technical Description	Max inlet pressure	10 bar	
	Pipe connection	External Diameter 15 mm	
	Pressure class	PN 16	
Specification, standard product	2 pcs gasoutlets for Medical Oxygen		
	 1 pc gasoutlet for Medical N₂O 1 pc gasoutlet for Medical Air 		
	 1 pc gasoutlet for Evacuation Ejector 		

Linde: Living healthcare

THE LINDE GROUP

AGA

Capacity and pressure drop


- The capacity for terminal units depends on the gas type, distribution pressure and accepted pressure drop.
 - Evacuation capacity: 25 l/min.

•

Gas type	Distribution Pressure kPa/bar	Max pressure drop kPa/bar	Capacity Ipm
Medical Oxygen	400/4	10/0,1	85
	500/5	10/0,1	85
Medical N ₂ O	300/3	10/0,1	80
	400/4	10/0,1	90
Medical Air	400/4	10/0,1	120
	600/6	10/0,1	130
Vacuum/Gas Evacuation	400/4	5/0,05	90
and Gas Outlet	400/4	5/0,05	90
Instrument Air	700/7	50/0,5	375
	800/8	50/0,5	400

ŝ

Measurements

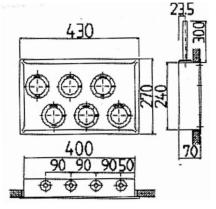

Recessed mour	nting:
Width	430 mm
Height	270 mm
Depth	70 mm

400		
	240	14 - 7

Exposed mounting:	
Width	400 mm
Height	240 mm
Depth	70 mm

Article Numbers	Art. No.	Denomination	
	301448	Recessed Gasoutlet Panel	
	301440	1440 Exposed Gasoutlet Panel	
	310918	Nozzle kit 0,4-0,7 (to change suction capacity)	

Contact us Linde Healthcare is committed to working with healthcare providers and regulatory authorities to continuously promote safe use of medical products and improve patient care.

We provide medical gas products, therapies, technical solutions and services to hospitals, clinics, nursing facilities, emergency management services and home healthcare providers around the world.

With our long experience and understanding of healthcare realities, you can depend on solutions that are delivered and serviced to the highest possible standards of quality, safety and efficiency.